



P: (732) 447-1100  
F: (732) 249-0008  
[www.pmdmfg.com](http://www.pmdmfg.com)


---

## LARYNGOSCOPES

---



*Expect the best in Laryngoscopy*

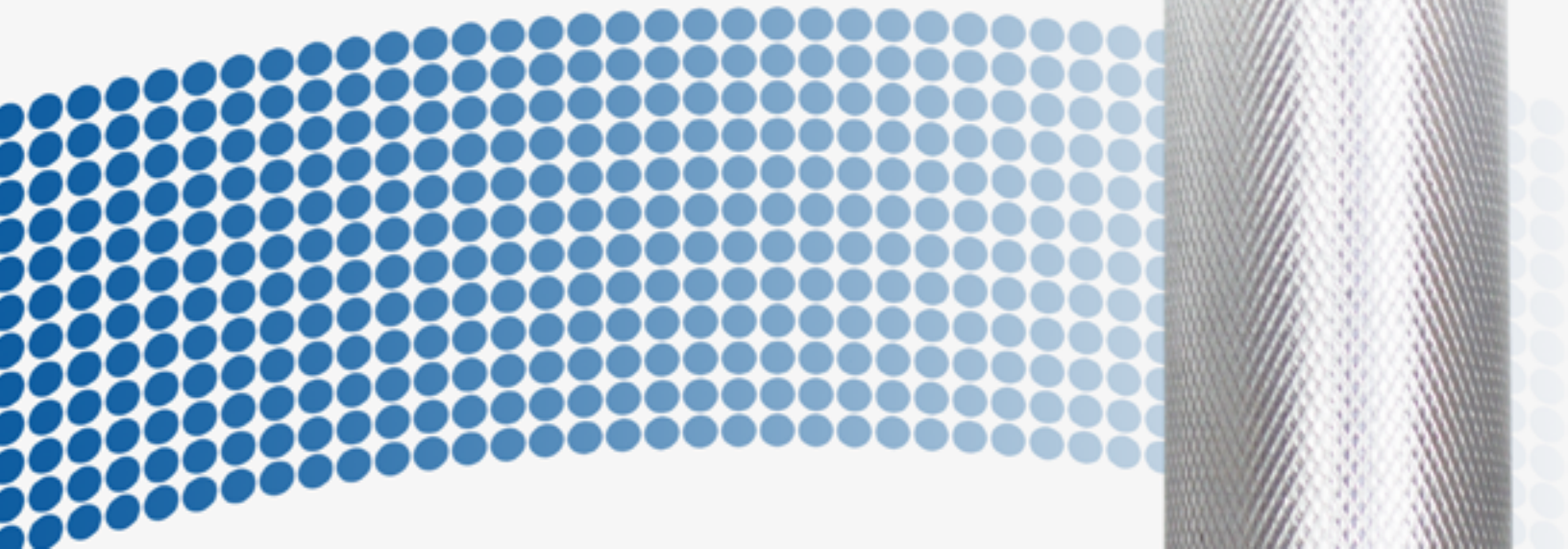

# ABOUT OUR PRODUCTS



## **PMD the laryngoscope specialists**

Whether you're looking for a traditional laryngoscope, or specialty device, we offer a tremendously wide portfolio of reusable and disposable blades.

We innovate to meet the clinician's needs. Here are some of the other products we offer.



# PMD

## WHAT SHOULD I KNOW ABOUT PMD?

Precision Medical Devices (PMD) is a leading manufacturer of high quality, great value laryngoscope blades and handles with a wide selection of the two most used blade shapes: Macintosh and Miller. Special sizes such as the Miller 00 are available. Our handles are manufactured of nickel coated brass with balanced heft and sure grip surface.

Our products use an LED lighting system providing strong, white light. Our Fiber Optic system is an excellent conduit of the illumination and this system is a traditional favorite of many anesthesiologists.

For those who prefer Conventional laryngoscopes, our fastest growing category, we offer a complete range of shapes and sizes.

Disposable or single-use blades have become the blade of choice for many institutions. They provide cost advantages, save re-processing time and expense, and help prevent cross contamination, a key objective of Joint Commission. We manufacture our disposable blades only from surgical grade stainless steel.

All Precision Medical laryngoscope blades and handles are manufactured under strict supervision to our exacting quality standards.



# LARYNGOSCOPE INTRODUCTION

## LARYNGOSCOPE BLADE SHAPES

We offer two standard blade shapes. Macintosh is the most widely used of the curved blade type and is the preferred shape for most adult intubations. Miller is straight, and designed to be longer and narrower than other blades. Most pediatric and neo-natal intubations are performed with Miller blades.

## FIBER OPTIC LIGHTING SYSTEMS

Fiber Optic transmits the illumination from our LED lamp in the handle through the fiber optic tube in the blade. The fiber optic tube is an excellent conduit of the light and this system is used by many anesthesiologists. Fiber Optic systems, handles, and blades are indicated by green marking.

## CONVENTIONAL LIGHTING SYSTEMS

Conventional lighting has the lamp in the blade. The light is activated when the contact point in the blade meets the contact point in the handle. PMD uses an LED lamp in its blades to provide a strong white light. The conventional system is growing rapidly in popularity.

## DISPOSABLE LARYNGOSCOPE BLADES

Disposable or single-use blades have become the blade of choice in many institutions. They provide cost advantages, save re-processing time and expense, and help prevent cross contamination. PMD has one of the broadest ranges of disposable surgical grade stainless steel blades available.

## EVERY BLADE

- Is made of AISI Surgical grade #304 Stainless Steel
- Fits International Standard (ISO 7376)
- Bright, yet cool focused LED Light

# SINGLE USE LARYNGOSCOPES

## A BRIGHT SOLUTION TO CROSS-CONTAMINATION

- LED Illumination Provides a whiter, brighter light for optimal illumination. Delivers 5X the light output of standard incandescent blades.
- 100% Stainless The ONLY Standard 100% Stainless Steel Single Use LED Blades and Single Use Laryngoscope Handles available.
- 100% Stainless steel blade smoothly releases from handle, with minimal effort.

## MINIMIZE RISK AND MAXIMIZE OUTCOMES

- Safety Addresses Joint Commission laryngoscope safety requirements. Single use blades / handles prevent cross-contamination /Improves reliability.
- Fixed light source - LED is fixed to blade for improved patient safety. LED is cool in operation – safer for neonates.
- Cost Effective Single use eliminates reprocessing/cleaning costs.
- LED utilizes minimal energy consumption, saving on battery replacement.
- Blades are compatible with standard conventional handles.





# FIBER OPTIC REUSABLE LARYNGOSCOPES

FIBER OPTIC REUSABLE LARYNGOSCOPES

## MACINTOSH PREMIUM BLADES



Size 0  
9311210



Size 2  
9311212



Size 1  
9311211



Size 3  
9311213



Size 4  
9311214

*Fits ISO Standard Handles*

Order Number	Size
9311214	4
9311213	3
9311212	2
9311211	1
9311210	0

# FIBER OPTIC REUSABLE LARYNGOSCOPIES

## MILLER PREMIUM BLADES



*Fits ISO Standard Handles*

Size 4  
9312214

Size 3  
9312213

Size 2  
9312212

Size 1  
9312211

Size 0  
9312210

Size 00  
9312200

Order Number	Size
9312214	4
9312213	3
9312212	2
9312211	1
9312210	0
9312200	00

# FIBER OPTIC DISPOSABLE LARYNGOSCOPES

FIBER OPTIC DISPOSABLE LARYNGOSCOPES

## MACINTOSH SINGLE-USE METAL BLADES



Size 0  
9311010



Size 2  
9311012



Size 1  
9311011



Size 3  
9311013



Size 4  
9311014

*Ask about our disposable handles*

Order Number	Size
9311014	4
9311013	3
9311012	2
9311011	1
9311010	0




# FIBER OPTIC DISPOSABLE LARYNGOSCOPIES

## MILLER SINGLE-USE METAL BLADES



*Ask about our disposable handles*

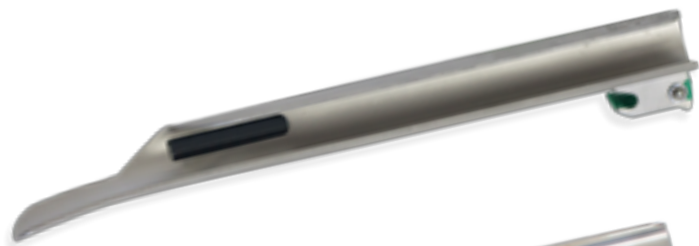

Size 4  
9312014



Size 3  
9312013



Size 2  
9312012



Size 1  
9312011



Size 0  
9312010



Size 00  
9312000

Order Number	Size
9312014	4
9312013	3
9312012	2
9312011	1
9312010	0
9312000	00

# CONVENTIONAL REUSABLE LARYNGOSCOPES

CONVENTIONAL REUSABLE LARYNGOSCOPES

## MACINTOSH PREMIUM BLADES



Size 0  
9301210



Size 2  
9301212



Size 1  
9301211



Size 3  
9301213



Size 4  
9301214



*Fits ISO Standard Handles*

Order Number	Size
9301214	4
9301213	3
9301212	2
9301211	1
9301210	0

# CONVENTIONAL REUSABLE LARYNGOSCOPES

## MILLER PREMIUM BLADES



*Fits ISO Standard Handles*

Size 4  
9302214

Size 3  
9302213

Size 2  
9302212

Size 1  
9302211

Size 0  
9302210

Size 00  
9302200

Order Number	Size
9302214	4
9302213	3
9302212	2
9302211	1
9302210	0
9302200	00

# CONVENTIONAL DISPOSABLE LARYNGOSCOPES

CONVENTIONAL DISPOSABLE LARYNGOSCOPES

## MACINTOSH SINGLE-USE METAL BLADES



Size 0  
9301010



Size 2  
9301012



Size 1  
9301011



Size 3  
9301013



Size 4  
9301014

*Ask about our disposable handles*



Order Number	Size
9301014	4
9301013	3
9301012	2
9301011	1
9301010	0

# CONVENTIONAL DISPOSABLE LARYNGOSCOPES

## MILLER SINGLE-USE METAL BLADES

*Ask about our disposable handles*



Size 4  
9302014

Size 3  
9302013

Size 2  
9302012

Size 1  
9302011

Size 0  
9302010

Size 00  
9302000

Order Number	Size
9302014	4
9302013	3
9302012	2
9302011	1
9302010	0
9302000	00



# LARYNGOSCOPE HANDLES

## LARYNGOSCOPE HANDLES

### CONVENTIONAL HANDLES

All PMD Handles are made according to ISO Standards

#### REUSABLE



Large  
9300003



Medium  
9300002



Small  
9300001



Stubby  
9300000

#### DISPOSABLE



One Size  
9306502

Order Number	Description	Size
9300003	Reusable	Large
9300002	Reusable	Medium
9300001	Reusable	Small
9300000	Reusable	Stubby
9306502	Disposable	One Size
9307582	Disp Metal	One Size

# LARYNGOSCOPE HANDLES

## FIBER OPTIC HANDLES

*PMD F/O handles always use an LED light source for a cool and bright white light*

### REUSABLE

### DISPOSABLE



Large  
9310003




Medium  
9310002



Small  
9310001



Stubby  
9310000

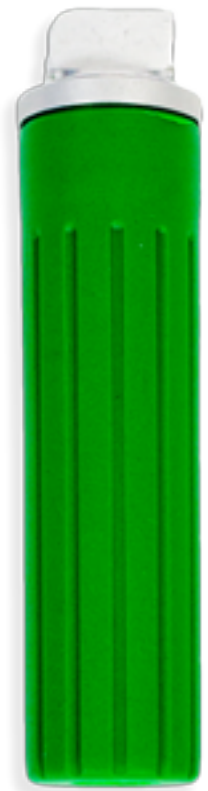

One Size  
9316502

Order Number	Description	Size
9310003	Reusable	Large
9310002	Reusable	Medium
9310001	Reusable	Small
9310000	Reusable	Stubby
9316502	Disposable	One Size
9307582	Disp F/O Metal	One Size



P: (732) 447-1100  
F: (732) 249-0008  
[www.pmdmfg.com](http://www.pmdmfg.com)

**Distributed By:**

© 2015 Cen-Med Enterprises, Inc. All rights reserved.